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SUN

Bluetooth-Körperfettwaage

Bluetooth Body Fat Scale

Bluetooth Vücut Analiz Baskülü

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Usage Instructions	26
Meaning of Symbols.....	26
Responsibility	26
CE Declaration of Conformity.....	26
Intended Use	26
Unauthorized Use.....	27
Safety Warnings	28
Opening the Box	29
User-friendly Tips	29
Features/Specifications.....	30
Performance Index.....	31
Battery Installation.....	31
Key Illustration.....	31
LCD Display.....	32
To Convert Weight Unit.....	32
Normal Weighing Mode (bypassing parameter setting)	33
Fat/Hydration/Muscle/Bone/Calorie Measuring Mode	34
Warning Indications.....	37
About Fat	38
About Calorie.....	38
How To Set Up And Use Our App.....	38
Advise For Use & Care	41
Faq	41
Fail to connect to your device.....	41
LCD has no responses	41
Advise for use & care	42
Device Storage.....	42
Recycling.....	42
Troubleshooting.....	43
Shipping	43

Dear Customer,

We hope you enjoy using your product. Thank you for choosing **Fakir Hausgeräte**, a leader in innovative technologies. Please read this user manual carefully before using the product and keep it for future reference.

Usage Instructions

These usage instructions are prepared to ensure that the device is used properly and safely. For your safety and the safety of others, please read the user manual before operating the device.

Meaning of Symbols

In this user manual, you will see the following symbols:

 **WARNING!**

This symbol indicates hazards that could lead to injury or death.

 **CAUTION!**

This symbol indicates hazards that could cause damage to your product or malfunction.

 **NOTE!**

Indicates important and useful information in the user manual.

Responsibility

Fakir Elektrikli Ev Aletleri Dış Tic. A.Ş. accepts no responsibility for any damage or repairs made outside Fakir Authorized Service Centers or for improper use of the device resulting from non-compliance with the instructions in this user manual.

CE Declaration of Conformity

 This device complies fully with the European Union's 2014/30/EU Electromagnetic Compatibility and 2014/35/EU Low Voltage Directives. The product bears the CE mark on the device type label.

Intended Use

The product model is designed for home use. It is not suitable for workplace or industrial use. Any malfunction resulting from such use will not be covered by our authorized service.

Unauthorized Use

In the following cases of unauthorized use, the device may be damaged or cause injuries:

- **Restricted Capabilities and Supervision:** This device should not be used by persons (including children) with limited physical, sensory, or mental capabilities or lack of experience and/or knowledge, unless they are supervised by a person responsible for their safety or given instructions regarding the safe use of the device. Children should be supervised to ensure they do not play with the device.
- **Use with Limited Abilities:** The device can be used by children aged 8 years and above and persons with limited physical, sensory, or mental capabilities or lack of experience and knowledge only if they are supervised, given instructions on safe use, and understand the potential hazards.
- **Packaging Materials:** Keep materials like bags included in the package away from children. There is a risk of suffocation!
- **Cleaning and Maintenance:** Children under 8 years should not clean or maintain the device unless they are supervised.
- **Device and Cable Safety:** Keep the device and its cable out of the reach of children under 8 years old. Children should be supervised to ensure they do not play with the device. Electrical devices are not toys for children. Therefore, use and store the device out of the reach of children. Children may not be aware of the dangers associated with the use of an electrical device. Ensure the power cord does not hang down from the device.

Safety Warnings

When using an electrical device, always follow the precautions below to prevent fire, electric shock, and personal injury:

- **Read the User Manual:** Carefully read the information in the user manual and keep it for future reference. If you give the device to a third party, be sure to include the manual.
- **First Use:** When you remove the product from its packaging, make sure all accessories and attachments are complete. If any parts are missing, contact your nearest Fakir Authorized Service Center. Keep the protective bags of the accessories and attachments away from children and infants to prevent the risk of suffocation.
- **Check the Condition of the Device:** If you suspect that there is any issue with the device, its casing, cable, or plug, do not operate the device and contact the nearest Fakir Authorized Service Center.
- **Cable and Plug Safety:** Do not use the device if the power cord or plug is faulty or damaged. Do not use the device with an extension cord and ensure that the power cord does not come into contact with heat sources or sharp edges. When unplugging, always pull the plug itself, not the cord. Do not pull on the cord or lift the device by the cord.
- **Child Safety:** Never leave the machine unattended while it is plugged in. Do not allow children to play with the device unsupervised. Store the device's power cord out of the reach of children. Electrical devices are not toys for children.
- **Use of the Device:** Do not use the device outdoors. Your machine may be adversely affected by dirt or rain. Do not immerse the machine in water and protect it from rain or moisture. The device operates with electrical power, which poses an electric shock hazard. Therefore, do not touch the plug with wet hands. Do not use the device to clean flammable or explosive materials (e.g., gasoline or alcohol), burning or smoldering objects (e.g., cigarettes, matches, or hot ashes).
- **Maintenance and Cleaning:** Before cleaning and maintenance, unplug the device and allow it to cool down.
- **Usage Restrictions:** Use the electrical device only as intended. Using the device contrary to the instructions may shorten its lifespan and pose safety risks.

CAUTION!

Carefully read the safety rules and usage instructions in this user manual.

Keep this user manual for future reference.

Opening the Box

- Before first use, remove the device and accessories from the box and check for damage or defects. If the device is damaged or defective, do not use it and contact Fakir Authorized Service.
- Repairs should only be performed by Fakir Authorized Service Centers.
- Repairs by unauthorized persons may be harmful to the user.

WARNING!

Do not use a faulty or defective product.

NOTE!

The values indicated on the device or in other printed documents provided with the product are obtained under laboratory conditions according to relevant standards. These values may change depending on the product's use and environmental conditions.

User-friendly Tips

1. During Measurement

- A. The device is designed to allow auto step on. Please always ignore the first reading and only start taking effective record from the second time weighing.
- B. Always remove your shoes and socks, wipe clean your feet before using.
- C. Step on the weighing platform gently.
- D. It is advisable to take measurement at the same time of the day.
- E. Measured result could be misleading after intensive exercise, excessive dieting or under extreme dehydration condition.
- F. Always weigh/use the scale on a hard and flat surface. When you are measuring, please keep still.
- G. The measured data of the following people may have deviation:
 - Children under 10 years old and adults over 100 years old. (only weight can be determined)

- Adult over 70.
- Body building athletes or other occupational athletes.

H. People with the following conditions should not use the device:

- Pregnant women
- Those with symptom of edema
- Those on dialysis treatment
- Those who use heart pacemaker or with other implanted medical device.

2. Safety Warning

- The Fat% Hydration% content evaluated by the Fat Monitor Scale is only for reference (not for medical purpose). If your Fat% Hydration% content is over/below normal level, please consult your doctor for more advice.
- Keep away from strong electro-magnetic field when using the scale.
- Don't use the scale on the damp floor, neither could stand on the edge to one side of the scale, otherwise it may tip!

Features/Specifications

1. Scale Body

- Using high precision strain gauge sensor
- Multifunction: test body fat %, body hydration %, body muscle%, body bone% and suggest calorie intake requirement.
- Built-in memory for 12 different users
- Optional normal weighing mode (bypassing parameter setting): Auto step on weighing and consecutive add-on weighings
- Low battery and overload indication

2. Bluetooth

- Directly check personal weight, body fat%, body hydration %, body muscle%, and body bone% information through graphs.
- Set password to protect personal information and store multi-user data;
- Simple and clear graphs current weight trends, form which you are able to keep track of the change;
- Bluetooth wireless technology;
- Connect with different smart devices respectively;

Performance Index

Index / Item	Weighing Capacity	Fat%	Hydration%	Muscle%	Bone%	Calorie
Index Range		4.0% ~60.0%	27.5% ~66.0%	20.0% ~56.0%	2.0% ~20.0%	
Division	0.1kg - 0.2lb	0.1%	0.1%	0.1%	0.1%	1 KCAL

Battery Installation

Name	Illustration	Installation
2x1.5V AAA Battery		AAA battery installation.
		Please install battery according to polarity.

NOTE!

Please use corresponding battery. Please remove battery if it is not in use in a long time (more than three months).

Key Illustration

SET (ON/SET) : Turn on the scale and prompt to next parameter. (Please refer to the Setting Personal Parameter overleaf)


▲: Value up by each interval press.

Press and hold the button to have consecutive value added.

▼: Value down by each interval press.

Press and hold the button to have consecutive value reduced.

LCD Display

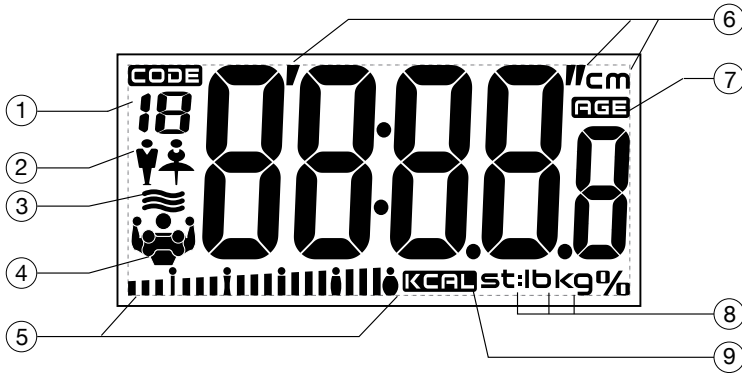- | | |
|--------------------------|-----------------|
| 1. User memory code | 6. Height Unit |
| 2. Gender | 7. Age |
| 3. Hydration% | 8. Weight Unit |
| 4. Fat%/Muscle%Bone% | 9. Calorie Unit |
| 5. Fat status indication | |


To Convert Weight Unit

If the scale is with both measuring system of Metric(kilograms) and Imperial (pounds, stones), you can convert weight unit as follow:

Condition 1

If the unit switch is on the back of the scale,


- Find unit switch on the back of scale.
- Pull up and down the unit switch to convert scale unit between kg-lb-st.

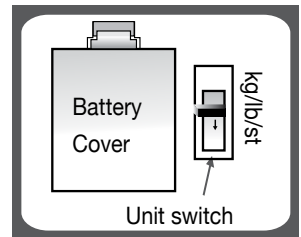Condition 2

If the unit switch is on the back of the scale,


- Find unit switch on the back of scale.
- Press the unit conversion button when scale is power on, LCD shows current weight unit. Press the button again to convert weight unit.

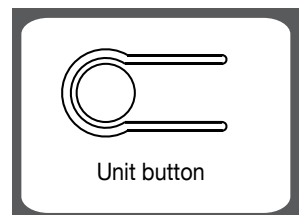Normal Weighing Mode (bypassing parameter setting)

1. Auto Step On Weighing

After installing battery, place the scale on flat floor. Keep stand still on scale when testing. Then wait for seconds, display shows the same weight reading for approximate a few seconds.

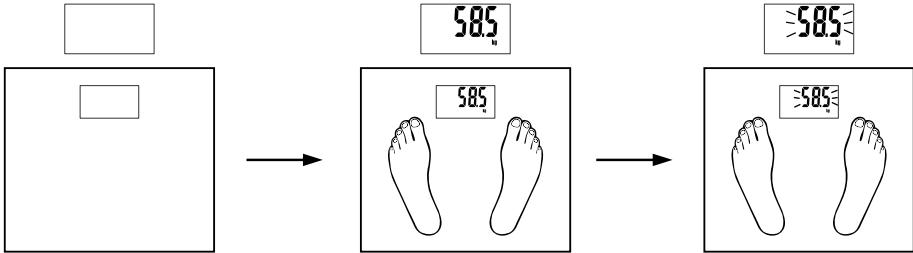Install battery and place the scale on flat floor. To make sure accuracy, place the scale on a hard, flat floor and avoid carpet and soft surface.


Step on scale slightly (Auto-on) for a while, wait for weight reading display (as shown above).

Reading result will flash after weight stabilized.

2. Consecutive Add-On Weighing

After digits stabilized, additional loadings can be conducted. Ensure each add-on loading has to exceed 2kg to get a new reading.

Example: If you want to get a weighing value of your baby, you can weigh yourself first, and then hold your baby together, without stepping off the scale. Then you will get the total weight reading. Subtracting the previous reading of yourself from the new reading, you will get the difference as your baby's weighing value.

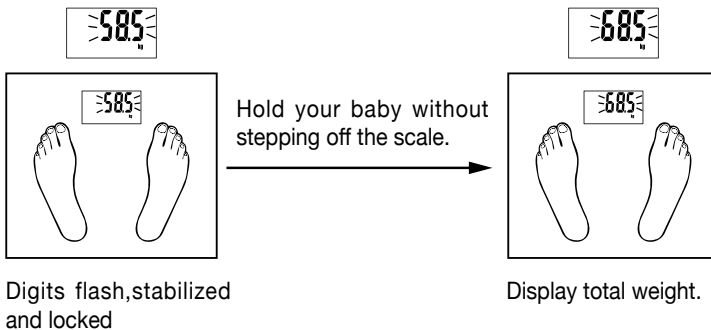Digits flash, stabilized and locked

Display total weight.

Fat/Hydration/Muscle/Bone/Calorie Measuring Mode

1. Setting personal parameters

At the first time when you install battery and switch on, all default parameter are set as follow:


Memory	Gender	Height	Age
1	Male	165cm/5'05.0"	25

Parameters can be set in the following range:

Memory	Gender	Height	Age
1-12	Male / Female	100-250cm 3'03.5"- 8'02.5"	10-100

Parameter Setting

Example: the 3rd user, female, 168cm, 28 years old, can set her parameters as follow:

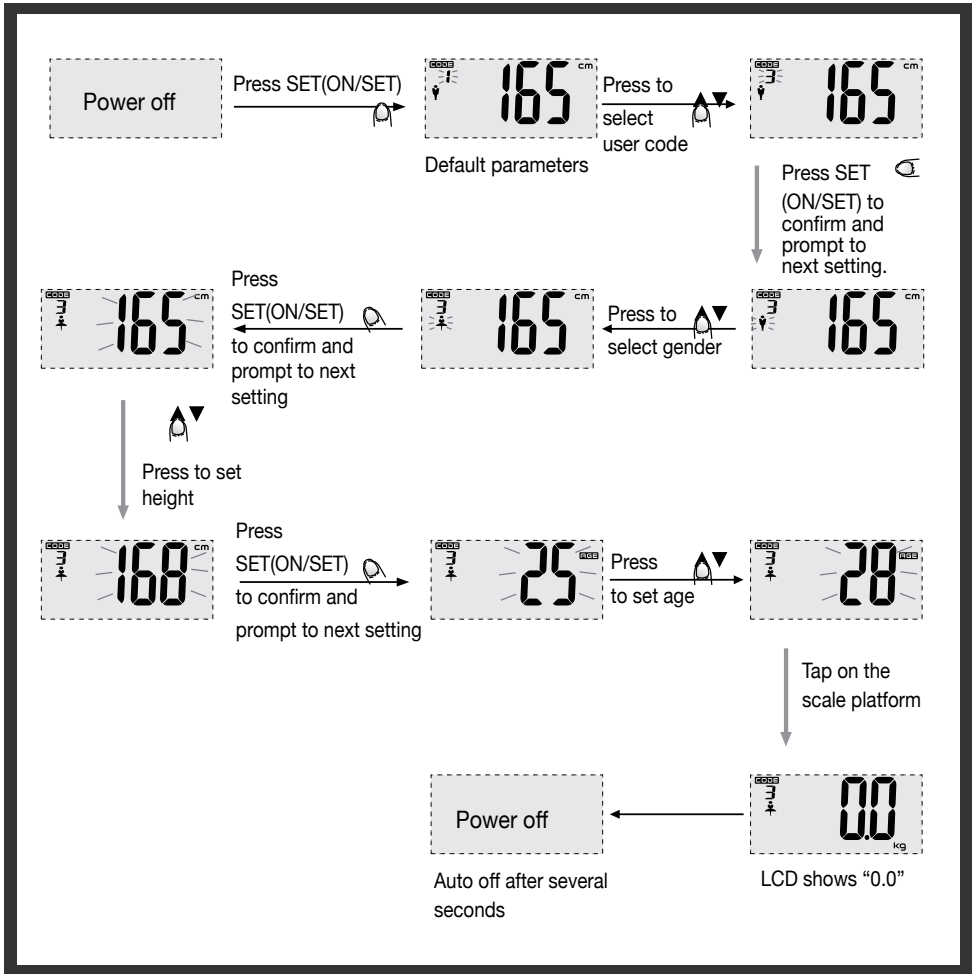 **NOTE!**

During parameter setting, if idle for 6 seconds, the parameters already set will be memorized automatically and LCD shows "0", meaning body fat measuring will be started.


 **NOTE!**

If LCD prompt to show "0" during parameter setting, press SET(ON/SET) for 3 seconds to return setting mode, then you can set user code, gender, height and age.

2. To Start Measuring


- (Make sure your feet are touching with metal electrodes on the scale platform, or body fat can not be measured.)
- To ensure accuracy, please remove your shoes & socks before stepping on.
- Example: the 3rd user, female, 168cm, 28 years old can measure as follow:

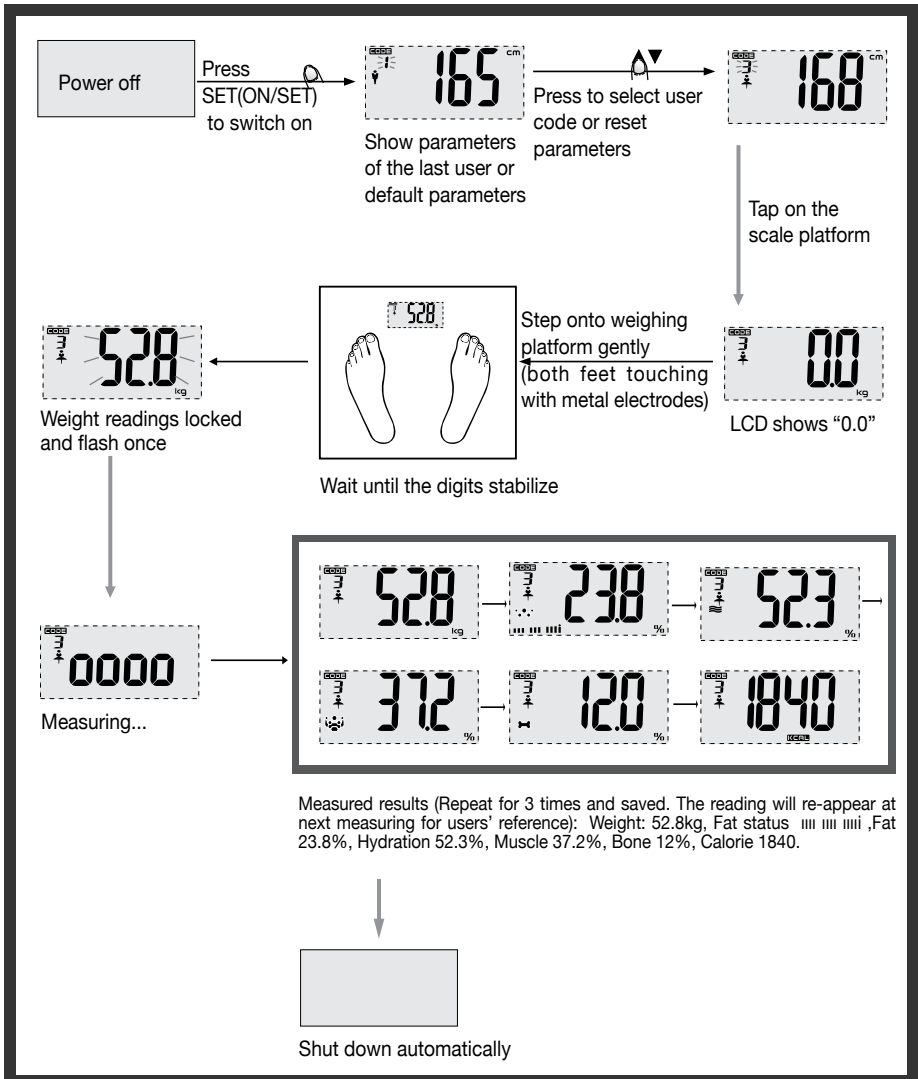Warning Indications

Low Battery Indication

The battery power is running low, please replace with a new battery.

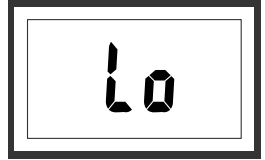Over-load Indication

The weighing subject on the platform exceeds the maximum capability of scale.

Please step off to avoid damage.

Low Fat% Indication

The Fat% is too low. Please lead a more nutritional diet and take good care of yourself.

High Fat% Indication

The Fat% is too high. Please watch your diet and do more exercise.

Re-test

Error exists, please re-test to get the correct result.

About Fat

1. About Fat%

- Fat% is an index % of fat content in human body.
- The fat content in human body is of a certain level, overabundance of fat will bring bad effects on metabolism, thereby, various illnesses may invade. Here goes an old saying, “Obesity is hotbed of illness”.

2. Working Principles

As electrical signal of certain frequency transmits through human body, the impedance of fat is stronger than that of muscle or other human tissues. Through the use of biological resistance analysis, this product will pass a safe frequency of low voltage electrical current through human body, thus, the Fat% will be reflected according to the resistance differences when the electrical signal transmitting through human body.

About Calorie

1. Calorie & Energy Consumption

Calorie: Energy Unit

The data reading indicates the energy you need to sustain the daily metabolism and ordinary activities required for your weight, height, age and sex.

2. Calorie & Losing Weight

Control the calorie intake combined with proper exercise is an effective way to lose weight. This is a very simple theory. As the daily energy intake can not support the energy consumption of the human body, the saccharide and fat stored inside will decompose and provide energy to human body.

That leads to weight loss.

How To Set Up And Use Our App

1. Get the Application

- Go to the “App Store” or “Google Play”, and search “MovingLife” to download our Fine-designed Application.

Start Here!

- You can also easily access our application by scanning the QR code located on the right-hand side.

2. Bluetooth Setting

- Make sure that the Bluetooth is already on.
- If not, go to the “Settings > Bluetooth” to turn it on.
- Tap on the platform of scale to turn it on.

NOTE!


Due to the new technology of Bluetooth 4.0, there is no need to pair the smart device with scale, it will be matched automatically within the Bluetooth signal coverage of 5-10 meters.

3. Registering and Connecting

Open APP MovingLife, connect the device with APP and submit your personal profile.(Before connection, please make sure Bluetooth is turned on and scale LCD is lighted up).

NOTE!

Please make sure your setted personal information is correct,in order to get accurate measurements.

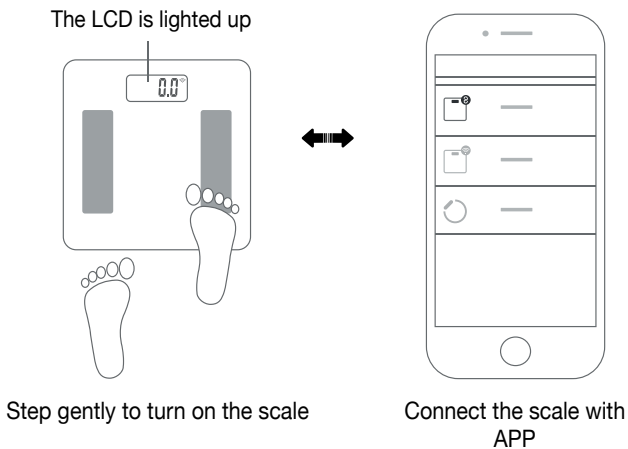4. Start Weighing

Place the scale on hard and flat floor. Stand on the scale with bare feet. The weight reading is stable and flashes,then begins to be transmitted. The scale shuts off after several seconds. APP will show and save all indexes and history record.

Place the scale on hard and flat floor.

Keep your bare feet on the electrodes of scale.

Check the measurement data on APP.

NOTE!

To make a proper measurement of body fat, please keep your bare feet on the electrodes of scale.

NOTE!

- The daily record for one user is always the last record measured and saved within day.
- The Fitness Assessment Chart is only for your reference, not determining everyone's real healthy status.
- There is also a Fitness Assessment Chart in the manual.
- You can modify your basic info, by clicking your avatar to the setting mode.

5. Tracking your weight and other indices

In order to track the trend of those indices, clicking the second icon on the bottom bar, you will be taken to the line graphs interface that shows the tendency of your weigh, fat%, hydration0/, Muscle%, Bone% and Kcal.

6. More Functions

- A. Click the third icon on the bottom bar.
- B. Customize languages, users management, weighing unit, or even setting a goal of weigh / bodyfat / BMI that it will make the tendency of the results be more referable.

ATTENTION!

The application is approving for your better experience, please keep your attention on the changes and improvements of the application with every update.

Advise For Use & Care

1. DO NOT use the scale in humid, overheat or over cold environment. (Temperature range from 5oC to 35oC)
2. DO NOT use the scale on soft mat, lawn, slippery floor or uneven floor.
3. DO NOT drop, shock or strike the scale as it is a precise instrument.
4. Keep the scale clean. DO NOT use acid cleaning agent.
5. DO NOT use the scale on a place with intense shake. Keep standing firmly on scale when in use.
6. To ensure accuracy, please test at the same time of the sequential days.

Faq

The application can't obtain weighing results from Bluetooth scale.

- Please confirm the Bluetooth icon on the application interface is in dark blue, showing connected successfully.
- Confirm the Bluetooth is on in your device.
- Confirm the Bluetooth scale is on within the valid range of signal coverage.
- Check the battery of Bluetooth scale, replace it if necessary.
- Restart your device and set up Bluetooth connection.

Fail to connect to your device

- Please make sure only one device is setting up connection with the Bluetooth scale.
- Please make sure your device is equipped with Bluetooth 4.0 or above.

LCD has no responses

- Please first check the batteries and replace them if necessary.
- Repairs may only be carried out by our authorized customer service or dealers.

Advise for use & care

1. The platform will be slippery when it is wet. Keep it dry!
2. Keep standing still during the course of weighing.
3. DO NOT strike, shock or drop the scale.
4. Clean the scale with a damp cloth and avoid water getting inside the scale. DO NOT use chemical/abrasive cleaning agents.
5. Keep the scale in cool place with dry air.
6. Always keep your scale in horizontal position.
7. If failed to turn on the scale, check whether battery is installed, or battery power is low. Insert or replace new battery.
8. If there is error on display or can not turn off the scale for a long time, please take off the battery for about 3 seconds, then install it again to remove soft-fault. If you can not resolve the problem, please contact your dealer for consultation.

Device Storage

- If you do not plan to use the appliance for a long while, store it in a cool and dry place.
- Ensure that the appliance is unplugged and its parts are completely dry.
- Ensure that the appliance is kept out of reach of children.

Recycling

When your product reaches the end of its life, follow these steps for proper recycling:

- Cut the device's cable to make it unusable.
- Dispose of the device in designated waste bins for such devices according to the applicable laws in your country.
- Electrical waste should not be disposed of with household waste.
- Throw the device into bins designated for electrical waste.
- The device's packaging is made from recyclable materials.
- Dispose of packaging materials in appropriate recycling bins.
- Dispose of the packaging materials in appropriate recycling bins.
- The device complies with the WEEE (Waste Electrical and Electronic Equipment) Directive.

NOTE!

- Check if you can fix the problem or defect yourself before calling customer service. For this, please refer to the “**Troubleshooting**” section. If you cannot resolve the issue or defect, call customer service.
- Before calling, please have the **model name** and product **serial number** ready. You can find this information on the **product's type label**.

Help Hotline: (+90 212) 444 0 241

Troubleshooting

Your device has undergone quality control to avoid any defects. However, if the device does not work for any reason, please try the solutions listed in the following table before contacting Fakir Authorized Service.

Shipping

To prevent damage to the device, follow these steps to transport your device safely:

- Transport the device in its original packaging.
- If the original packaging is not available, use a suitable and soft packaging that protects the device.

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